


AP Connect Editions

DOWNLOAD
FULL 60-DAY DEMO

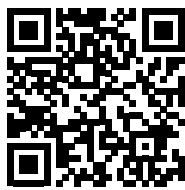

www.anton-paar.com/
apc-demo

	Start (free)	Standard	Pharma
	↓	↓	↓
	For all customers who want to familiarize themselves and grow with AP Connect, and want to connect a single instrument.	Broad-based. For labs and applications across industries.	Focused on labs in intensely regulated environments. For operators with specific stringent testing requirements.
Especially for these industries	<div><div>- Food & beverages</div><div>- Petroleum</div><div>- Chemical</div><div>- Flavors & fragrances</div></div>	<div><div>- Food & beverages</div><div>- Petroleum</div><div>- Chemical</div><div>- Flavors & fragrances</div></div>	<div><div>- Pharmaceutical</div><div>- Personal care</div><div>- Healthcare</div><div>- Nutraceutical</div><div>- Government compliance</div></div>
License	Free	One-time license fee or annual subscription fee	One-time license fee or annual subscription fee
Instrument connectivity	1 instrument	Up to 30 instruments	Up to 30 instruments
Data collection	✓	✓	✓
Data analysis (results browser)	✓	✓	✓
Data export	✓	✓	✓
Configurable automatic exports	1 rule	Up to 100 rules	Up to 100 rules
Supported file formats for data export	PDF, CSV	PDF, XLS, CSV	PDF, XLS, CSV
Desktop installation	✓	✓	✓
Server installation	✓	✓	✓
Client installation (multi-user access)	1 client	✓	✓
User management supporting Active Directory	✓	✓	✓
Fileless data transfer (REST API) ¹		✓	✓
Connect to instruments from other vendors ²		✓	✓
Control of instruments (automated measurement series)		✓	✓
Roles and permissions			✓
Review & approval process			✓
Electronic signature			✓
Audit trail			✓
Compliant with Title 21 CFR Part 11			✓

As per the features shown in the table, compliance with the ALCOA++ data integrity principle is ensured.
¹ Via separate license (extra license costs apply).
² Additional adapters needed.

© 2023 Anton Paar GmbH | All rights reserved.
Specifications subject to change without notice.
D96IP008EN-D

Software for Your Paperless Lab

AP Connect



AP Connect.

Let the Data Flow.

Data can be complex, but this is easy.

Connect and collect, plug and play. All your lab data – just a click away. AP Connect is a lab execution software that connects instruments from Anton Paar and other vendors, communicates measurement information, and ensures compliance. It's your data hub in the lab – your world right at your fingertips.



Paperless.

Data at Your Fingertips.

Why is AP Connect different than anything else? The versatile system architecture allows a plug-and-play integration within a few minutes. No complex, time-consuming implementation. No bells and whistles. Just time-saving functionality and increased productivity.

Tens of thousands of measurements in one digital space. 100 % efficiency. Fun, even. Catalyze your lab work.

Reward your data – in three easy steps

- 1

INSTALL
Install AP Connect and, optionally, the instrument adapters in as little as four minutes
- 2

LOG IN
Log in to AP Connect in seconds with your Windows access credentials (Active Directory)
- 3

REGISTER
Register your instruments in less than one minute



←
Do my instruments work with AP Connect?
See all of the 60+ Anton Paar instruments that connect.
See all supported instruments from other vendors
connected via instrument adapters

www.anton-paar.com/apc

Connect.	Collect.	Check.
<div>→ Plug and play</div> <div>→ Connect instruments from Anton Paar and other vendors</div>	<div>→ One central database for all your measurement data</div> <div>→ Access from across the company network</div> <div>→ Work with information instead of searching</div>	<div>→ Call up your data anytime to make faster decisions</div> <div>→ Audit trail function for full traceability in your lab work</div> <div>→ Configurable approval process of all your data with electronic signature</div>

Zero transcription errors ensure data quality.

Zero paperwork saves time and frees up resources for data analysis.

All data for analysis in one place enables fast and safe decisions.

Control.	Communicate.	Comply.
<div>→ User management using the Active Directory</div> <div>→ Secure assignment of roles & permissions</div> <div>→ Trigger measurements with 10+ parameters directly from your PC</div>	<div>→ One single interface for all your instruments to your data management system</div> <div>→ Remote data access: for fileless data transfer via the standardized REST interface</div> <div>→ File-based automatic and manual exports</div>	<div>→ The combination of all AP Connect features ensures full compliance in regulated environments (e.g., FDA 21 CFR Part 11).</div>

Batch measurement increases sample throughput.

Configure your communication paths once, and your data will flow freely.

System audit trail and permission management guarantee full traceability.